


The *secure* solution for Trendelenburg Positioning

Simple. Secure. Reusable.

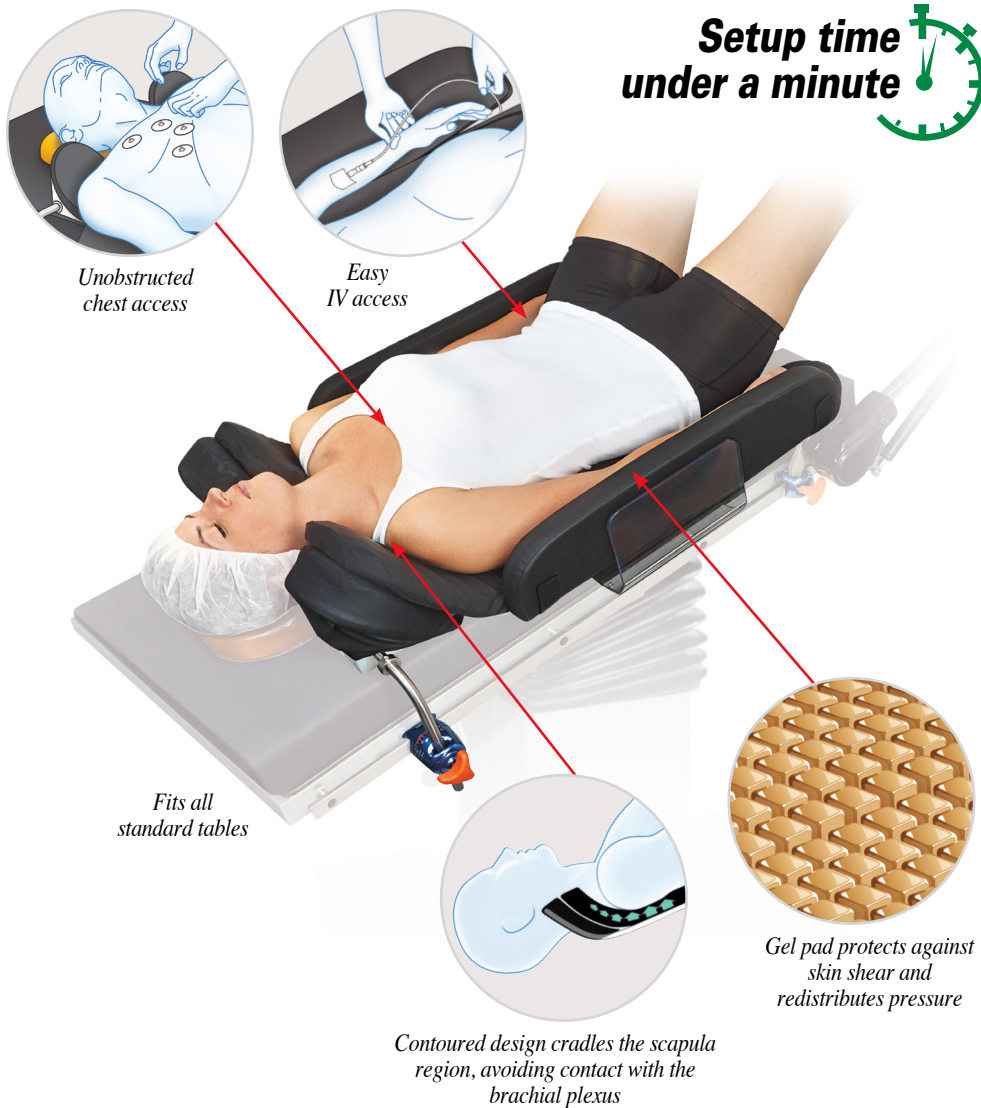

**Setup time
under a minute**

■ Cradles patient to maximize pressure redistribution

■ Unobstructed chest and IV access

■ Improved design for easier airway access

SecureFit[®] TPS

Secure, easy-to-use solution dramatically reduces setup time, resists slipping, evenly redistributes pressure, and is easy to clean.

Stabilizing Wedges with improved design are easier to remove for faster airway access (see back for details).

Engineered to minimize slipping Cradling the body's natural curvature, the SecureFit TPS reduces peak focal pressure over the shoulder and minimizes slipping with no ancillary shoulder braces required.

Maximum pressure redistribution The **TPS Gel Pad** conforms to the patient's anatomy and redistributes pressure to provide support.

IV access simplified Access IV lines quickly and easily through convenient openings between the arm support and gel pad.

See reverse side for more details.




Specifications:

TPS Foundation: 20"W x 17"L x 6"H • **TPS Gel Pad:** 22½"W x 39"L x 1½" thick (1" gel surface with ½" Polyfoam base) • **Perineal cutout:** 9"W x 13"L • **TPS Stabilizing Wedge Pads:** 1" memory foam over Polyfoam • **Cover:** Black anti-shear elastic • **Antimicrobial** • **Fluid-proof fabric** • **Fluid-resistant sewn seams** • **CAL 117-compliant** • **Stabilizing Posts:** 9¼"L x ¾" diameter • **Stainless steel** • **Toboggan:** 10"L x 8"H x 10¼"D • **Toboggan Arm Supports:** 25"L x 4½"H Polyfoam pads

TPS Arm Supports and Toboggans

Integrated, easy-to-clean vinyl-covered foam pads work with toboggans to protect ulnar nerve while providing easy IV access to the arms.

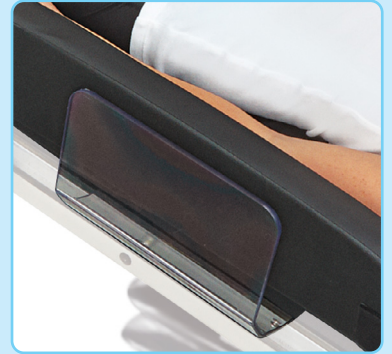

Ordering Information:

#2970015137 SecureFit TPS Trendelenburg Positioning System

Replacement Items:

- #2970015138** TPS Gel Pad
- #2970015139** TPS Stabilizing Wedge Pads, pair
- #2970015135** Shoulder Wedges, pair
- #2970015140** Stabilizing Posts, pair
- #2970015141** Toboggans, pair

MD • 65 California residents, see below.

Patent #9,283,133



TPS Foundation

Strong, durable, nonskid foundation is contoured to cradle patient's natural curvature, and locks onto the surgical table with clamps for a secure hold.

Stabilizing Wedges

Improved design! Wedges slide straight out for easier removal and faster access to the airway.



Watch the video and learn more at [AliMed.com/securefit](https://www.ali-med.com/securefit)

Or scan here:



Any order placed before 3 p.m. EST, if in stock, ships the same day!
Prices, product availability, and specs subject to change.

California Proposition 65 Statement

65 **WARNING:** This product can expose you to chemicals including Di(2-Ethylhexyl) Phthalate (DEHP), which is known to the State of California to cause cancer and birth defects or other reproductive harm. For more information go to <https://www.p65warnings.ca.gov>.



Pub. #7253REV1223

©2023 AliMed® Inc. All rights reserved.