


A necessity, not a luxury, for larger patients

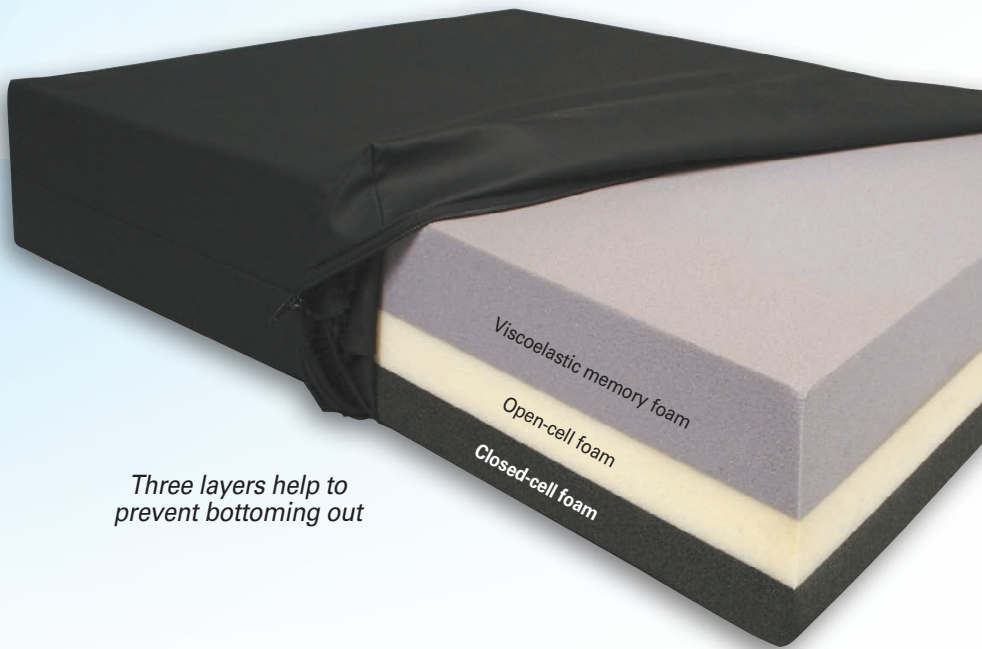

Three layers help to
prevent bottoming out

Best for patients who:

- ✓ Use wheelchairs up to 30"W
- ✓ Weigh up to 1,000 lbs.
- ✓ Are less mobile and at increased risk for pressure injuries



Need help choosing the right wheelchair cushion?
Read our blog at [AliMed.com/choosing-a-cushion](https://www.ali-med.com/blog/choosing-a-cushion)

■ Three layers of foam designed for bariatric users' needs

■ 2"H top layer of viscoelastic memory foam provides exceptional pressure relief

■ Base layer of closed-cell foam prevents bottoming out


Bariatric Wheelchair Cushions

Specialty seat cushions for bariatric wheelchair users are a necessity, not an accessory. Bariatric wheelchair cushions aid in distribution and reduction of pressure, help prevent bottoming out, and promote better posture and a straighter spine.

Unlike standard cushions that are made bigger but not stronger for bariatric patients, **Bariatric Tri-Foam Wheelchair Cushions** have three distinct layers of support—a 2"H top layer of viscoelastic memory foam for exceptional pressure relief, a 1"H middle layer of open-cell foam for pressure reduction, and a 1"H base layer of sturdy closed-cell foam to help reduce pressure and prevent bottoming out. Available in sizes up to 30"W to accommodate both weight and size, fulfilling your largest users' needs.

See reverse for
ordering information 

Bariatric cushions for maximum support

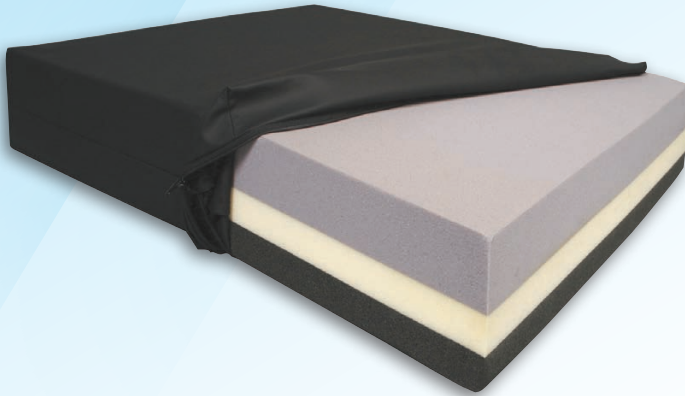

Bariatric Tri-Foam Cushions Layered foam fits wheelchairs up to 30"W without bottoming out.

4"H • Removable, washable cover • 1,000-lb. capacity

ITEM	W x D	ITEM	W x D
#2970000080	20" x 18"	#2970000085	26" x 20"
#2970000081	22" x 18"	#2970000086	28" x 18"
#2970000090	22" x 22"	#2970000087	28" x 20"
#2970000082	24" x 18"	#2970000088	30" x 18"
#2970000083	24" x 20"	#2970000089	30" x 20"
#2970000084	26" x 18"		

MD • 65 California residents, see below.

High-density pressure relief cushions

Medial-lateral thigh positioners for optimal posture

Removable fluid-proof stretch cover included with both models

Elite model has moderate front-to-back slope

High-density foam base


Shear-relieving gel at critical pressure points

Elite Cushions Multiple layers of gel and foam for optimal pressure reduction. **AT (Anti-Thrust)** model includes pre-ischial bar for greater pelvic stability.

5"H tapers in rear • Meets CAL 117 • 2-year warranty • Includes removable, launderable, black, fluid-proof stretch cover • 350-lb. capacity

ITEM	W x D
ELITE	
#2970000234	16" x 16"
#2970000235	18" x 16"
#2970000236	18" x 18"
#2970000237	20" x 18"
#2970000238	22" x 18"
#2970000239	24" x 18"
ELITE AT (ANTI-THRUST)	
#2970000240	18" x 16"

MD • 65 California residents, see below.

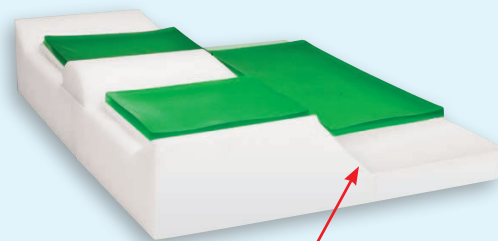

Anti-Thrust (AT) model with pre-ischial bar helps prevent forward sliding

Visit [AliMed.com/cushions](https://www.ali-med.com/cushions) to shop our full line.

Any order placed before 3 p.m. EST, if in stock, ships the same day!
Prices, product availability, and specs are subject to change.

California Proposition 65 Statement

65 ⚠️ **WARNING:** This product can expose you to chemicals including Di(2-Ethylhexyl) Phthalate (DEHP), which is known to the State of California to cause cancer and birth defects or other reproductive harm. For more information go to <https://www.p65warnings.ca.gov>.

Pub. #7672REV0324 Product #11524

©2024 AliMed® Inc. All rights reserved. Printed in USA.